

ZOLL AED 3[®] Trainer

Training AED and Accessories



Complete ZOLL AED 3 simulation with Real CPR Help real-time CPR feedback

The ZOLL AED 3[®] Trainer is the only automated external defibrillator (AED) training solution that provides Real CPR Help[®], real-time feedback on a trainee's compression quality. The trainer will allow students to learn proper rate and depth of CPR chest compression in a training setting.

The ZOLL AED 3 Trainer is designed to closely mimic the operation of the clinical AED, offering students a more realistic rescue experience. The ZOLL AED 3 Trainer includes:

- An active LCD display with full-color rescue images that mirror what a student will see on a clinical ZOLL AED 3 device
- Contact sensors in the training pads that detect pad contact with the manikin, eliminating the need for a remote control to support rescue simulations
- An on-device pause button that allows the instructor to pause a training scenario
- Long-lasting consumables, which save time and money
- Easy USB software update/language configuration

2015 AHA and ERC Guidelines Ready

Easy configuration: All configuration and training scenario management is performed on the trainer.

Configurable options include:

- Training scenario
- Operating mode (Semi-automatic or fully automatic)
- Lay-rescuer prompts
- Check-breathing prompts
- Voice-prompt volume
- Language (varies by part number ordered)

The ZOLL AED 3 Trainer includes everything an instructor needs:

- ZOLL AED 3 Trainer device
- ZOLL AED 3 Trainer CPR Uni-padz[®] electrode training harness (minimum of 100 training sessions)
- Trainer electrode gels/liner for ZOLL AED 3 Trainer CPR Uni-padz (minimum of 50 training sessions)
- Four (4) D cell batteries
- Operator's Guide

ZOLL[®]



Real-time CPR feedback with a large 2-color compression depth indicator



Configure device operation and rescue scenarios directly on the device

ZOLL AED 3 Trainer Specifications

Size: (D x W x H) 5.0 in x 9.3 in x 9.7 in
(12.7 cm x 23.6 cm x 24.7 cm)

Weight: 3 lbs 15oz. (1.8 kg) with batteries

Power: Four (4) D cells (included) for Training AED (Good for more than 16 hours of training)


Rescue Scenarios

Instructors can set up various rescue scenarios directly on the trainer using soft keys.




ZOLL AED 3 Trainer Accessories

Replacements must be ordered separately.

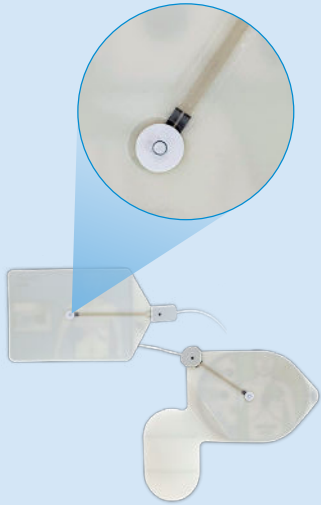

Pad sensors detect contact with the manikin, eliminating the need for a separate remote

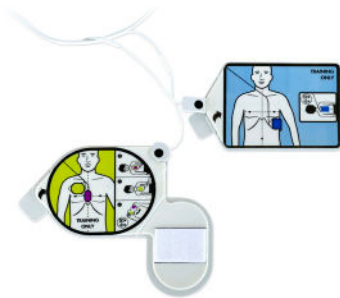

Trainer CPR Uni-padz[®] electrode replacement adhesive gels
Part# 8028-000012 (case of 5)



Trainer CPR Uni-padz electrode replacement liners Part# 8028-000013



Trainer demo manikin
Part# 8000-0835-01



Trainer USB Software Upgrade Kit
Part # 8028-000011-XX



CPR Uni-padz electrode training harness
Part# 8028-000010

ZOLL Medical Corporation
Worldwide Headquarters
269 Mill Road
Chelmsford, MA 01824
978-421-9655
800-804-4356

For subsidiary addresses and fax numbers, as well as other global locations, please go to www.zoll.com/contacts.

Specifications subject to change without notice.

©2019 ZOLL Medical Corporation. All rights reserved. Real CPR Help, Uni-padz, ZOLL, and ZOLL AED 3, are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners. The ZOLL AED 3 defibrillator is not available for sale in the United States, Canada, or Singapore.

Printed in U.S.A. MCN PP 1902 0365

ZOLL